


HOW TO USE YOUR *HEM HANS* FREQUENCY DEVICE


SWITCHES USED TO CONTROL THIS DEVICE.

- I. UP
- II. DOWN
- III. BACK
- IV. ENTER
- V. FAVOURITE
- VI. DELETE FAV
- VII. START/STOP
- VIII. PAUSE/CONTINUE

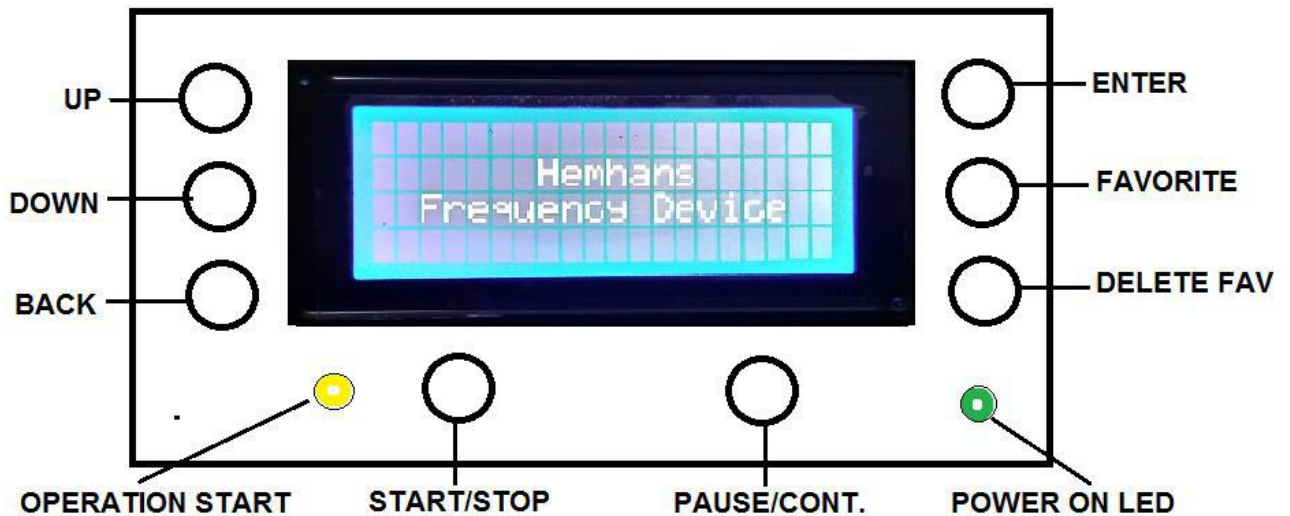

DEVICE OPERATION MANUAL

1. There is red on/off switch given on the right side to initialize the machine.
2. When we start the machine message "HemHans Frequency Device" will appear in the display.

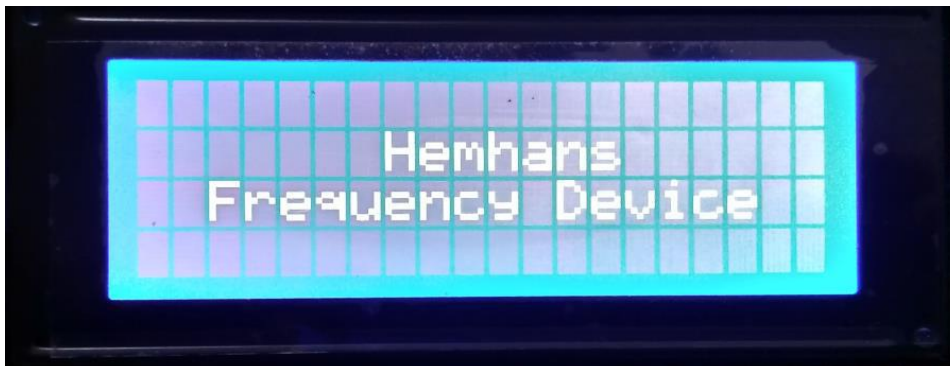

3. Press ENTER. There will appear two categories Old Freq & New Freq. To select one of the categories press UP or DOWN button.



4. After choosing the freq, press ENTER then you can see the diseases are classified in A to Z classes. Select one of the class presses UP or DOWN button. To get back to previous function press "Back" button.



5. Press ENTER on one of A - Z, we can access diseases of that alphabet. Select particular disease presses UP or DOWN button. Suppose if we want to take frequency for "Acne", He will have to enter in "A" class where all frequencies for Acne are available. Display shows "Disease


name - Acne” and “Sequence number [A001]”.To get back to previous function press “Back” button.



6. Press “START/STOP” button to start the frequencies. Now display shows disease name - “Acne”, Initial frequency - “F: 0.05” ,Total time in minute (T) - “T”:30”, Remaining time in minute (R) - “R”:29.At the same time ON/OFF LED will blink. The frequency will be changed automatically after 3 minutes. When we press “START/STOP” button again the operation will be stopped and device will stop to generate frequencies. At the same time ON/OFF LED will OFF.




7. Press “Pause/Cont.” button to pause the operation at any time. When we press this button display shows “Pause!!” and device will stop to generate frequencies. At the same time ON/OFF LED will continuously in ON state. When we press Pause/Cont.” again the device will resume its operation.

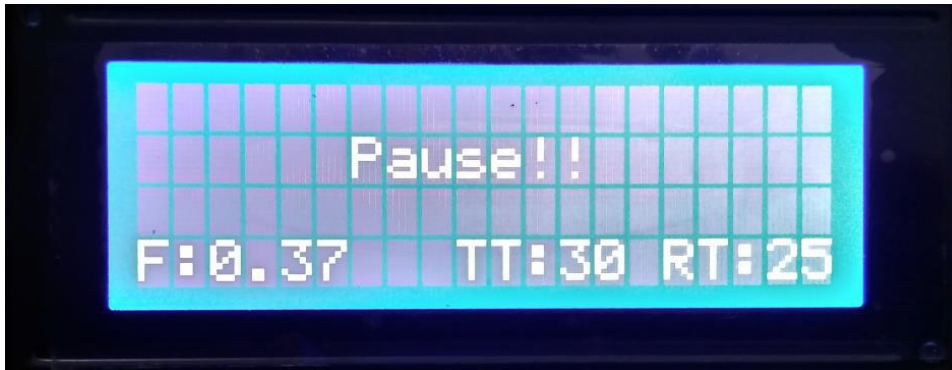


8. FAVOURITE: This switch is used to set your routine frequencies which you want to take frequently.



FAVOURITE: Suppose you want to take 'ACNE' frequently - Navigate to 'Acne' from the A-Z class and then click on Favourite button. When you click on Favourite button you will see the below list: -

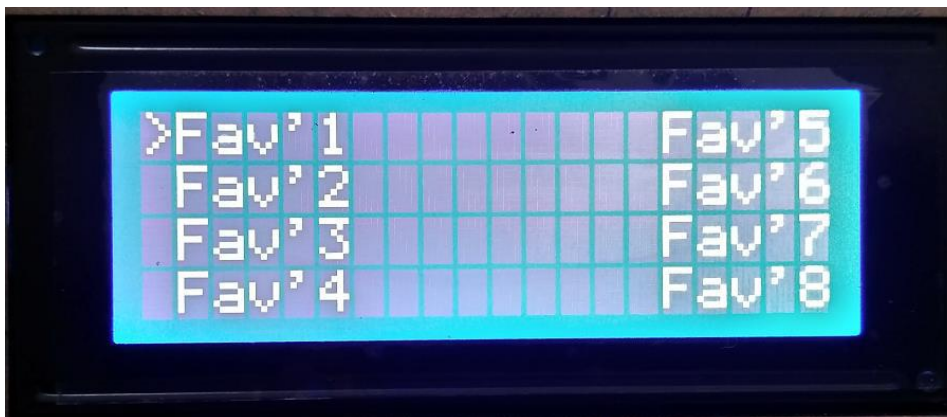

You can enter up to 8 favourite lists. Select any list using UP/DOWN button and click ENTER. You will see the below screen.

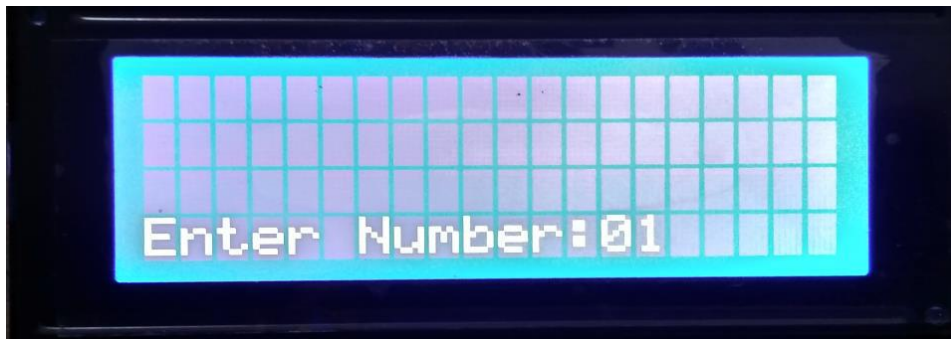

In each list you can enter 10 favourite diseases. Using UP/DOWN key, select a number and click ENTER. If you want to keep Acne fifth in the list, Enter Number: 05 and then click ENTER.

You will get a message 'Added Successfully' Now if you want to take the diseases directly click on FAVOURITE button and select the Fav'1 and START taking them.

9. DELETE + UP: You can skip a frequency by pressing DELETE and UP button together.
10. BATTERY: You can see four LED lights in the BATTERY section. When only one LED is lighted please charge your machine.
11. MUSIC: There is a Music Selection switch (3-Way Black Switch) at the right side of the machine, by using this selector switch you can select any one category (General mantra or Jain Mantra) out of two categories of Mantra. If you want to stop the music put the same switch to middle position. There is a Volume roller at the right side of the machine. By rotating it you can change the volume of the music. Using the red switch on the right side of the machine, you can change the mantra's that are played.

WARNING



Do not use Frequency Device if you are wearing a pacemaker or if you are pregnant.

Do not over use or abuse the Frequency Device.

Do not use while in water.

Always seek the advice of a medical professional with regards to any health issue.

DISCLAIMER

HEM HANS Frequency Device is based on the innovation made by Dr. Royal Rife and Dr. Hulda Clark both. The frequencies of bacteria/viruses are given by Dr. Royal Rife and by Dr. Hulda Clark in her book “The Cure for All Diseases”. Their Research is based on the observations and case-studies of the humans. We do not give any guarantee for the results and it may vary person-to-person.

The information given in this Handbook is for educational, informational, and investigational purposes only and is not to be construed as diagnosis of disease, treatment of disease, prevention of disease, or as a replacement for consulting a qualified health care practitioner. Be careful when investigating this technology! Protocols for this technology may need to be modified, or this technology may be contraindicated entirely, if you have a heart condition, are wearing a pacemaker, are pregnant, are nursing, have blood clots, are taking strong medications such as chemo, are wearing metal implants or stents, have breast implants, are especially sensitive to radio frequency (RF) or other electromagnetic radiation, have problems with your immune response, or have especially sluggish eliminative functions (colon, kidneys, liver and lymph system). Before using any equipment, and to see if you should even be experimenting with this technology, please read about these circumstances, and the precautions to take. The author, publisher and distributor are not responsible or liable for the results of your experimenting with rife technology or using any of the other therapeutic modalities described in this book. The reader accepts full responsibility for any and all consequences of his or her experiments with rife technology or any other electro medicine modality. If you have a medical condition, see a qualified health professional of your choice.

WARRANTY:

HEM HANS Frequency Device is warranted for 12 months from the date of purchase for manufacturing defects only. We do not undertake any warranty for the external or internal damages caused for mishandling of the same.

For any issues please contact: -

Dr. Hemant Broker
'Lal Bungalow',
Opposite Anjali Jewellers,
C.G. Road,
Ahmedabad - 380006
Gujarat, India

Phone No: - 07926401444, 07926405434